

# SCREENING FOR PESTICIDE RESIDUES IN CANNABIS USING A HIGH RESOLUTION ACCURATE MASS GC/Q-TOF

Philip L. Myle, Mei Wang, 2 Mahmoud A. ElSohly<sup>2,3</sup>  
 Ikhlas Khan<sup>2</sup> Chandrani Majumdar<sup>4</sup> and Mohamed Radwan<sup>4</sup>  
<sup>1</sup>Agilent Technologies, Wilmington DE; <sup>2</sup>National Center for Natural Products Research, <sup>3</sup>Department of Pharmaceutics and Drug Delivery, School of Pharmacy, University of Mississippi, <sup>4</sup>EISohly Laboratory Inc.

NACRW  
2018  
Poster # 67



## Introduction

**Cannabis is:**  
 Considered to be a schedule 1 drug by the US government same as heroin, LSD, ecstasy...  
 Legal for medical and recreational use in 9 states & DC  
 Legal for medical use in 29 states  
 Regulations are state-by-state  
 No uniformity in pesticide regulations

California legalized recreational use – Jan. 1, 2018  
 Canada will legalize recreational use this summer

California is the state with the most comprehensive testing requirements

Black market competes with legalized cannabis pesticides widely used on black market cannabis



## Results and Discussion

**Cannabis Extracts are Dirty! High Concentrations of Terpenes and Cannabinoids**

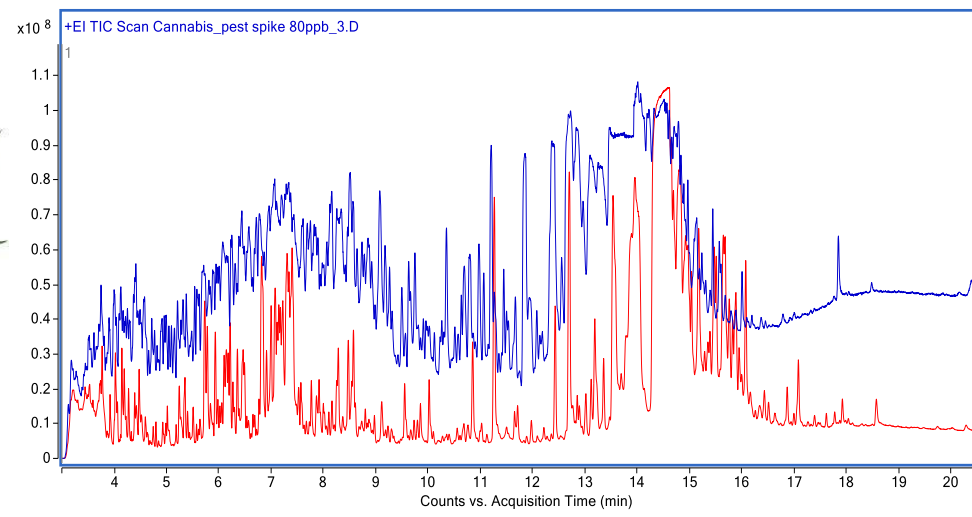


Figure 4. Chromatograms of two different cannabis extracts. Red: used procedure shown in Figure 2. Blue: used the modifications described below Figure 2.

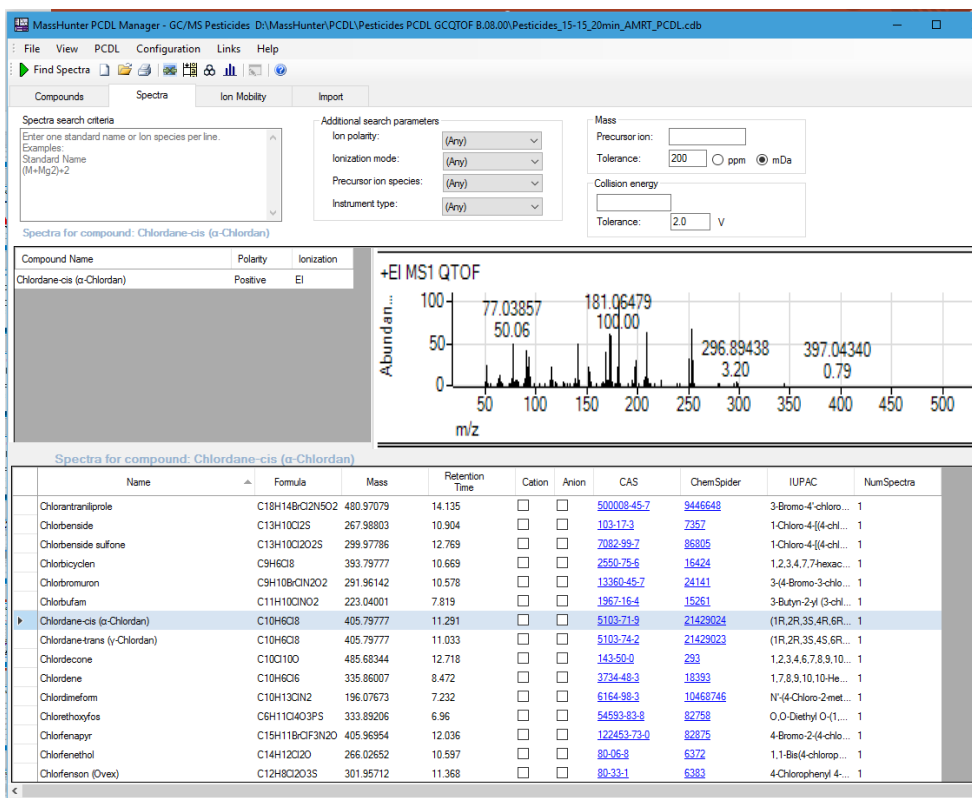


Figure 5. Exact mass Pesticide Personal Compound Database and Library containing HRAM spectra for >850 pesticides and related contaminants.


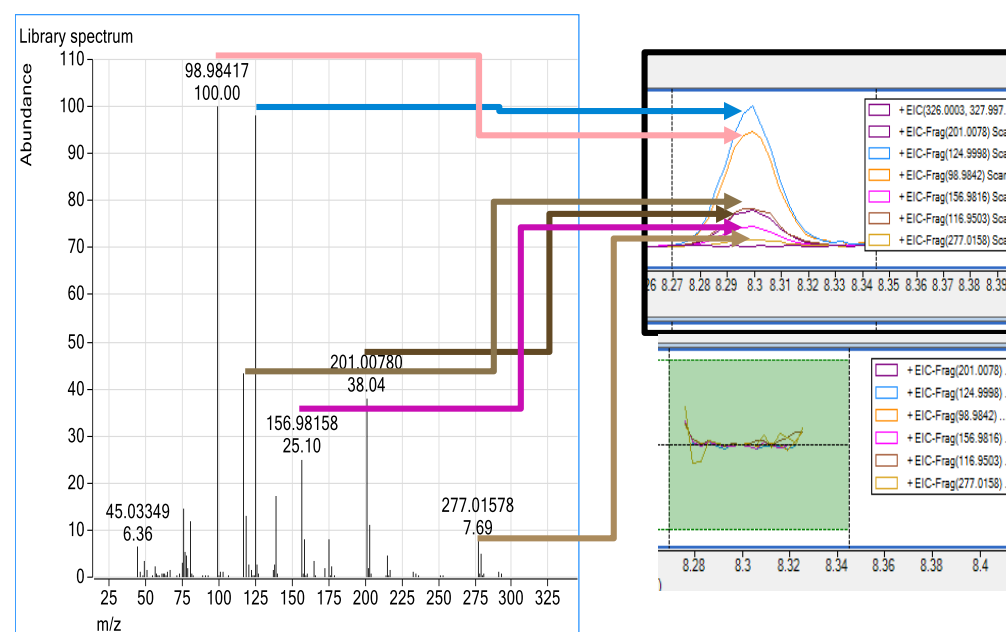


Figure 6. Software extracts six exact mass ions at the locked RT for each compound in the PCDL. EICs with the same RT and peak shape are "qualified." Other information is available: RT difference from PCDL, # of qualified ions, fragment ratio score, mass accuracy, and more.

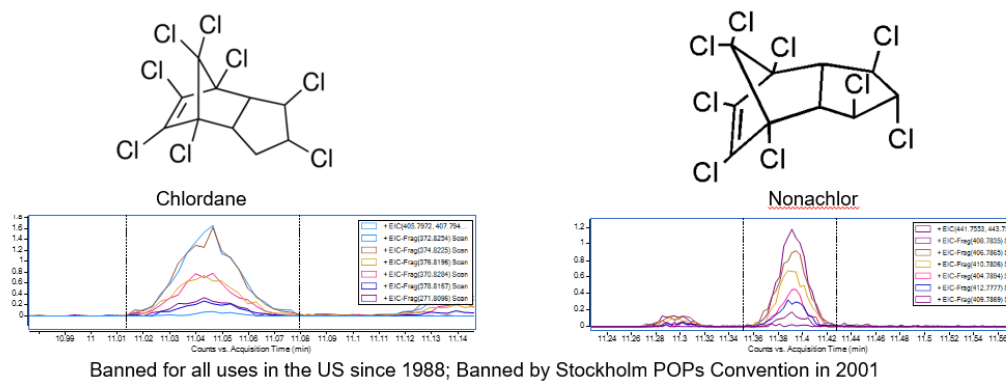

Banned for all uses in the US since 1988; Banned by Stockholm POPs Convention in 2001  
 Figure 7. Banned pesticides found in two confiscated samples

**Pesticides & contaminants found in 16 confiscated cannabis samples**

Pesticides & Metabolites	Pesticides & Metabolites
DEET	Benzyl Benzoate
Cypermethrins	gamma-Cyhalothrin
Tributyl phosphate*	2-Phenylphenol
Acrinathrin	Cyphenothrin
Jasmone	Pentachloroisole (Pentachlorophenol metabolite)
Empenthrin	
Dichlorvos	<b>Fire Retardants</b>
trans-Chlordane	Tris(3-chloropropyl)phosphate*
cis-Chlordane	Tris(beta-chloropropyl)phosphate*
trans-Nonachlor	
Diphenylamine	
2,6-Diisopropynaphthalene*	<b>PAHs</b>
Terbucarb	Phenanthrene*
Fenpropathrin	Fluoranthene*
Mylobutanil	Pyrene*
Carbaryl	<b>Yellow highlight = on California's Forbidden List</b>
1-Naphthol (Carbaryl metabolite)	<b>Blue: California has a Max Residue Limit</b>

## Quantitative Analysis

**Spiking samples**

Cannabis grown without pesticides at U. Mississippi Spiked with ~200 pesticides @ 10, 20, 40, 60, 100 & 200 ng/g before extraction

Chromatogram of cannabis extract shown in Figure 4 (blue)

Matrix interferences interfered with many pesticides

Reasonable calibration curves for 83 pesticides

R<sup>2</sup> for most >0.99

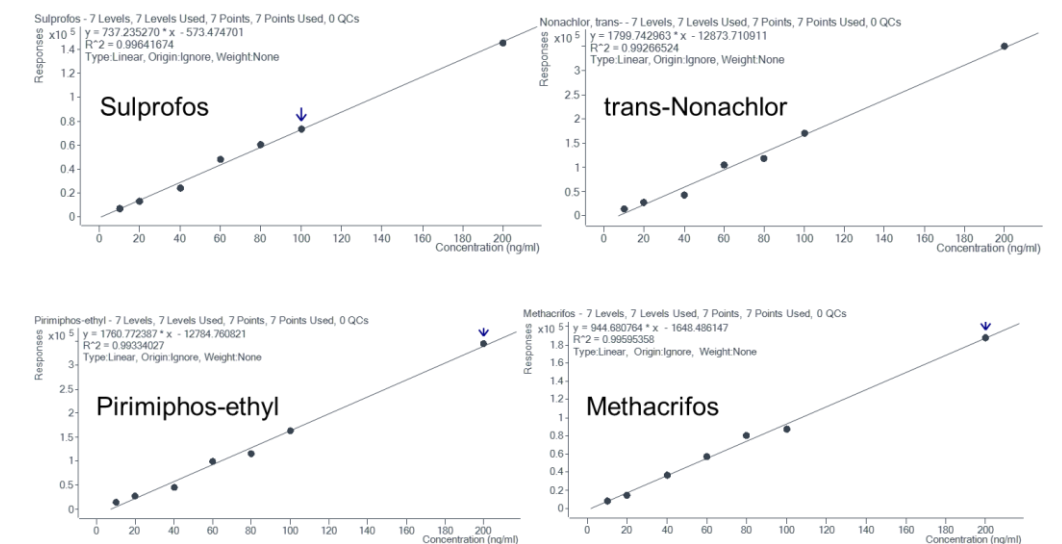


Figure 8. Calibration curves for select pesticides analyzed by GC/Q-TOF. Interferences from coextracted endogenous compounds prevented many pesticides from being quantified.

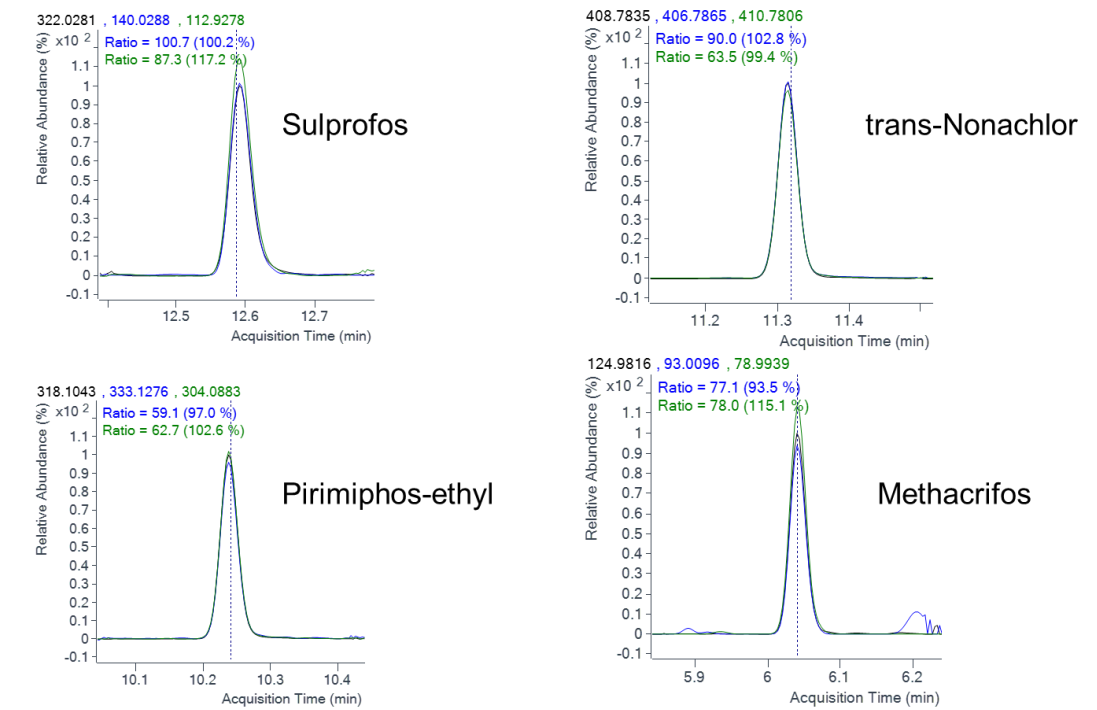


Figure 9. Quant and qualifier ions for selected pesticides analyzed by GC/Q-TOF

## The Role of GC/QQQ

Cannabis testing labs use GC/QQQ for pesticide residue testing because of its selectivity and lower cost.[1] One can quantify more than 200 pesticides but one only "sees" those compounds on the target list

## Next Steps for Cannabis Analysis by GC/Q-TOF

- Improve extraction process
- Dilute extracts by at least 100X
- Analyze samples using an Agilent 7250 GC/Q-TOF with
  - Improved sensitivity
  - Much higher resolution
  - Larger linear dynamic range
  - Low energy ionization
- Try MS/MS with the Q-TOF

## Conclusions

### Cannabis Extraction

- Two methods tried
  - QuEChERS extraction with 10:1 dilution
  - Modified QuEChERS with 3:1 dilution
- Both methods resulted in an extract with a lot of matrix
- Will try other extraction strategies with 100:1 dilution

### Qualitative Analysis

- 22 different pesticides identified in 16 confiscated marijuana samples
- OP fire retardants and PAHs are ubiquitous in the environment and may be environmental or processing contamination
- Some pesticides were probably missed because of matrix

### Quantitative Analysis

- Method developed for 82 pesticides
- Heavy matrix interfered with many pesticides
- Better cleanup and sample dilution will probably help

## References

[1] L. Asanuma, D. Miller, R. Jordan, M. Churley, and A. Macherone, A novel comprehensive strategy for residual pesticide analysis in cannabis flower, Agilent Application Note 5991-9030EN.

## Experimental

### Extraction Procedure: Modified QuEChERS



Figure 2. Extraction procedure

A second set of extractions used 1 g of cannabis & 3 mL each of acetonitrile & water & an additional back extraction with a second aliquot of water. This gave a much dirtier extract (Figure 4).

GC Parameters	
Injection volume	1 µL
MMI Temperature	280°C (Splitless)
Columns	Two 15 m X 0.25 mm x 0.25 µm HP-5MS UI
He Carrier gas	Col 1 = 0.979 mL/min; Col 2 = 1.179 mL/min
Oven Temp Program	60°C (1 min), 40°C/min to 170°C (0 min), 10°C/min to 310°C (3 min)
RT Locking	Chlorpyrifos locked to 9.954 min
Backflushing	3 min post run at 310°C and -7.87mL/min
Q-TOF Parameters	
Mode	TOF
Acquisition range	40 – 550 u
Acquisition rate	5 Hz
Collision gas	N <sub>2</sub> at 1.5 mL/min

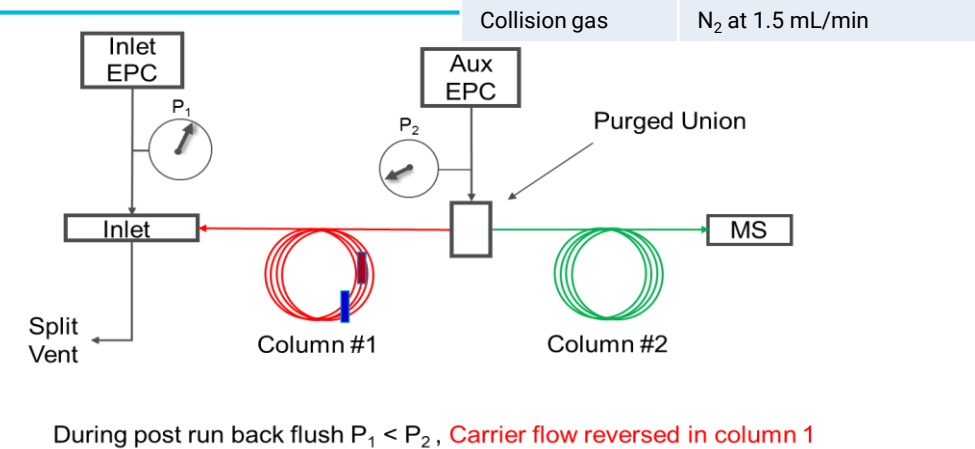


Figure 3. GC configuration for backflushing column 1