


Film Holder Accessory

for the Cary 4000, 5000, 6000i, 7000
UV-Vis-NIR spectrophotometers

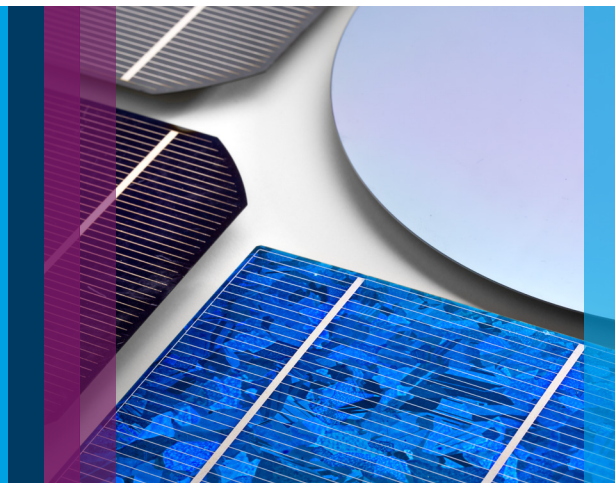


Automated multi-filter and film homogeneity measurements

The film holder accessory is designed to be mounted on the sample transport accessory and can be used to test the homogeneity of a sample, measuring the %T or Abs as a function of distance. The accessory provides a vertical mounting capability for films, gels, sheets of material or multiple filters.

The accessory is equipped with a wide variety of mounting options, with a large number of mounting holes providing flexibility in sample size. Mylar strips along the accessory protect delicate samples from scratching. The extended sample compartment accessory is required (sold separately).

Specifications

Maximum sample size	160 x 160 mm
Minimum distance interval	1 mm over 160 mm of travel.

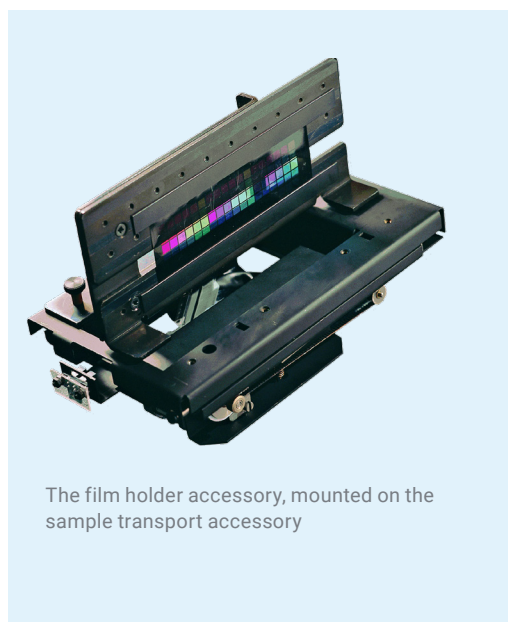

The film holder accessory, mounted on the sample transport accessory

For more information.

www.agilent.com/chem

DE44462.6551851852

This information is subject to change without notice.

© Agilent Technologies, Inc. 2021
Published in the USA, September 27, 2021
5991-1718EN