


GC CAPILLARY COLUMN APPLICATIONS: FLAVORS/FRAGRANCES

Free Fatty Acids (in Methanol)

COLUMN: 007-FFAP, Bonded
FFAP
15M. x 0.32mm I.D. x
0.25 μ m film
Cat. No.: 007-FFAP-15W-0.25F
Temperature: 140°C
Injector: 220°C
Detector: 250°C , FID
Carrier Gas: 25 cm/sec., Helium

1. acetic acid
2. proprionic acid
3. butyric acid
4. isovaleric acid
5. hexanoic acid

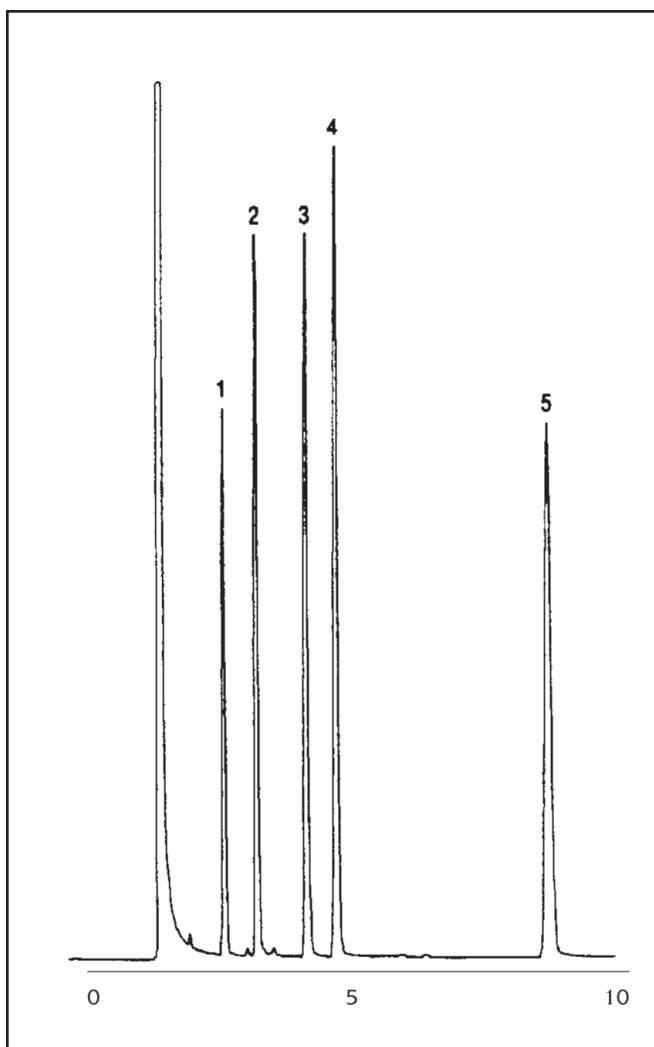

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com